

Oxygen ULR CHEMets® Kit

K-7511/R-7511: 0 - 20 ppb

Sampling

The most critical part of any dissolved oxygen test is sampling. The sample stream must be completely leak-free. To accomplish this, the sampling tube is vertically mounted with a tube of inert material connecting the sample point to the bottom of the sampling tube. Use stainless steel, type 304 or 316, or glass tubing with short neoprene connections. Do not use copper tubing, long sections of neoprene or other polymeric tubing.

Test Procedure

1. To remove trapped air bubbles, the system should be purged with water that is flowing at the fastest possible rate, and has a temperature of 180 - 210°F (80 - 100°C). New sampling systems should be purged for several hours, while those used routinely may require only a few minutes. **When the system is fully purged, reduce the flow to 500 - 1000 mL per minute and cool the sample to ambient temperature.**

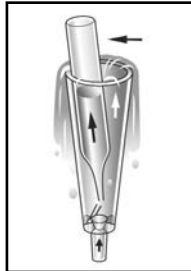


Figure 1

2. Place the ULR CHEMet ampoule, tip first, into the sampling tube. Snap the tip. The ampoule will fill leaving a bubble for mixing (fig. 1).
3. Gently invert the ampoule, allowing the bubble to travel from end to end **10 times**. Dry the ampoule.
NOTE: The bubble will travel most rapidly if the ampoule is held at a 45° angle while mixing.

4. **Immediately** obtain a test result by gently placing the ampoule, flat end first, into the comparator. Hold the comparator up toward a source of light and view from the bottom. Rotate the comparator until the best color match is found (fig. 2).

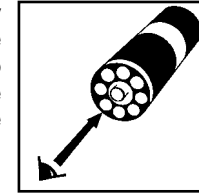


Figure 2

Test Method

The Oxygen ULR CHEMets®¹ test kit employs the Rhodazine D™ Method.^{2,3,4,5} Dissolved oxygen reacts with the pale yellow colored leuco form of Rhodazine D to produce a deep rose color. The resulting color is proportional to the dissolved oxygen concentration in the sample.

1. ULR CHEMets is a registered trademark of CHEMetrics, Inc. U.S. Patent No. 3,634,038
2. Rhodazine D methodology was developed by and is a trademark of CHEMetrics, Inc.
3. ASTM D 5543 - 15, Low Level Dissolved Oxygen in Water
4. ASTM Power Plant Manual, 1st ed., p. 169 (1984)
5. Department of the Navy, Final Report of NAVSECPHILADIV Project A-1598; Evaluation of CHEMetrics Feedwater Dissolved Oxygen Test Kit (1975)

Safety Information

Read SDS (available at www.chemetrics.com) before performing this test procedure. Wear safety glasses and protective gloves.

Visit www.chemetrics.com to view product demonstration videos.
Always follow the test procedure above to perform a test.



Simplicity in Water Analysis

www.chemetrics.com
4295 Catlett Road, Midland, VA 22728 U.S.A.
Phone: (800) 356-3072; Fax: (540) 788-4856
E-Mail: orders@chemetrics.com

Feb. 18, Rev. 15